



## Mini Pin Micro-Needle Protocol



Main Unit

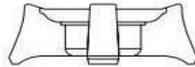

Stand

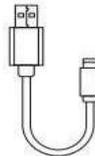

Charger




Cartridge Tip & Cap

### How to use:


1. Cleansethe skin thoroughly. Wewrecommend using **CLINICARE Concentrated Cleansing Foam**. Mix a small amount of the foam with water and make a light foam. Massage the skin and remove with lukewarm water and pat dry. Repeat if necessary.
2. Use a new and unused cartridge for every treatment. With clean hands, unseal the cartridge from its blister pack and separate the cap from the cartridge tip by turning counter-clockwise.
3. Fit just the tip cartridge into the main unit by aligning the notches and place it on the stand.
4. Fill the cartridge with a suitable mesotherapy solution. Wewrecommend **CLINICARE LMWHASerum**. Do not add more than 3ml at a time.
5. Securely fasten the cartridge cap onto the cartridge tip by rotating it clockwise.
6. Turn the Mini Pin on and check the battery. Press and hold grey power button for 2 seconds - if light blinks it will need to be re-charged for 30 minutes. Press power button on again briefly to switch the device on.
7. Check the flow of the serum before using by pressing the white button in the middle of the needles in cartridge tip 8-10 times.
8. Roll the Mini Pin slowly and smoothly over the face (or scalp). Do not apply pressure to the skin. Do not treat same area twice. Ensure device is at right angles to the skin.
9. Press grey power button twice to turn the Mini Pin off.
10. Apply a calming and soothing sheet mask such as **CLINICARE EGF Mask** for about 20 minutes. Massage in the rest of the mask.
11. Apply **CLINICARE SunShield Cream SPF30** or **Dermo Corrective Cream SPF50**. Sunscreens should be applied every day for at least one week after use.



Cleanse the skin thoroughly

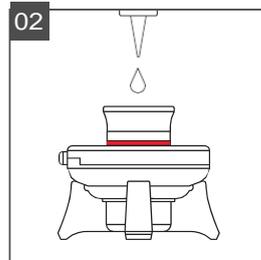

Place Mini Pin on stand & pour mesotherapy serum into the Cartridge Tip

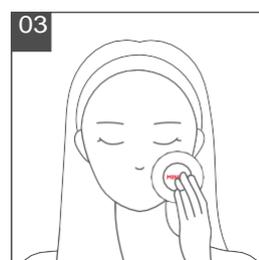

Roll the Mini Pin slowly over skin



Apply a soothing sheet mask

**Precautions:**

1. Please read instructions and precautions thoroughly before using the Mini Pin.
2. The cartridges are single use only and should be safely disposed of immediately after use. We are not responsible for any product or customer issues caused by reuse.
3. Ensure you use the recommended brands or hypoallergenic products suitable for at-home micro needling.
4. Avoid sun exposure on the day of use. A sunscreen with an SPF30 or above should be applied daily for at least one week after use.
5. Refrain from extreme temperatures such as saunas, sunbeds etc for at least 1 week.
6. Do not swim in chlorinated water for 24 hours.
7. Avoid electrolysis, waxing, bleaching (face), depilatory creams, laser hair removal on the area for at least 72 hours.
8. Do not use more than once a week.
9. Avoid using around the lips and eyelids.
10. Keep the product out of reach of children and pets.
11. Slight redness, tenderness and dryness may occur after use. These reactions are temporary and typically resolve within 3-4 days as the skin returns to normal. Refrain from rubbing or picking the skin.
12. If using it while charging, charge the Mini Pin directly from the outlet (5V/1A). Connecting using a USB port from a PC or Notebook is prohibited.

**Use is prohibited in the following cases:**

1. Metal allergy / Infectious skin disease (such as cold sores) / Uncontrolled diabetes / History of keloid scarring / Skin cancers / Haemophilia / Inflammatory acne, rosacea, eczema or psoriasis within the application area
2. Pregnant women / Breastfeeding.
3. Consult us prior to use if you have had recent ablative dermal procedures or aesthetic treatments.
4. If you notice any abnormalities such as skin rashes, inflammation, or allergies during use, please discontinue use and consult your GP.